

DISPOSABLE LINEAR STAPLER WITH KNIFE



QUALITY CARE THROUGH INNOVATION



Product Code	Description	Staple Line Length	Cutting Line Length	Staple Height	Staple Formation	Pack Qty
CAAB-40	Disposable Linear Cutter Stapler	43mm	40mm	4.8mm	2.0mm	3/box
CAZB-40D	Reloading Unit for Linear Cutter Stapler	43mm	40mm	4.8mm	2.0mm	1/box

Linear Stapler with Knife

LOW ANTERIOR RESECTION

Because Accessibility Matters!



Because accessibility matters!



Suitable for lower colorectal surgery: it combines horizontal stapling and vertical cutting and supports multiple firing, making operations more flexible.



Panther Healthcare Medical Equipment Co., Ltd
Beijing | New York | São Paulo | Toronto | Mumbai

E-mail: orders@pantherhealthcarecanada.com
Tel: 1-833-293-8971
45 Midpark Crescent, London, ON N6N 1A9
pantherhealthcarecanada.com



www.pantherhealthcarecanada.com

LINEAR STAPLER WITH KNIFE

*The smallest head and the largest jaw aperture on the market.
Facilitates access and mobility to place the stapler and the tissue exactly where it's needed*

Longer shaft makes the device easier to reach bottom of pelvis to conduct low rectal surgery

Straight cut edge and staple line - Well-formed staple line ensuring good blood circulation

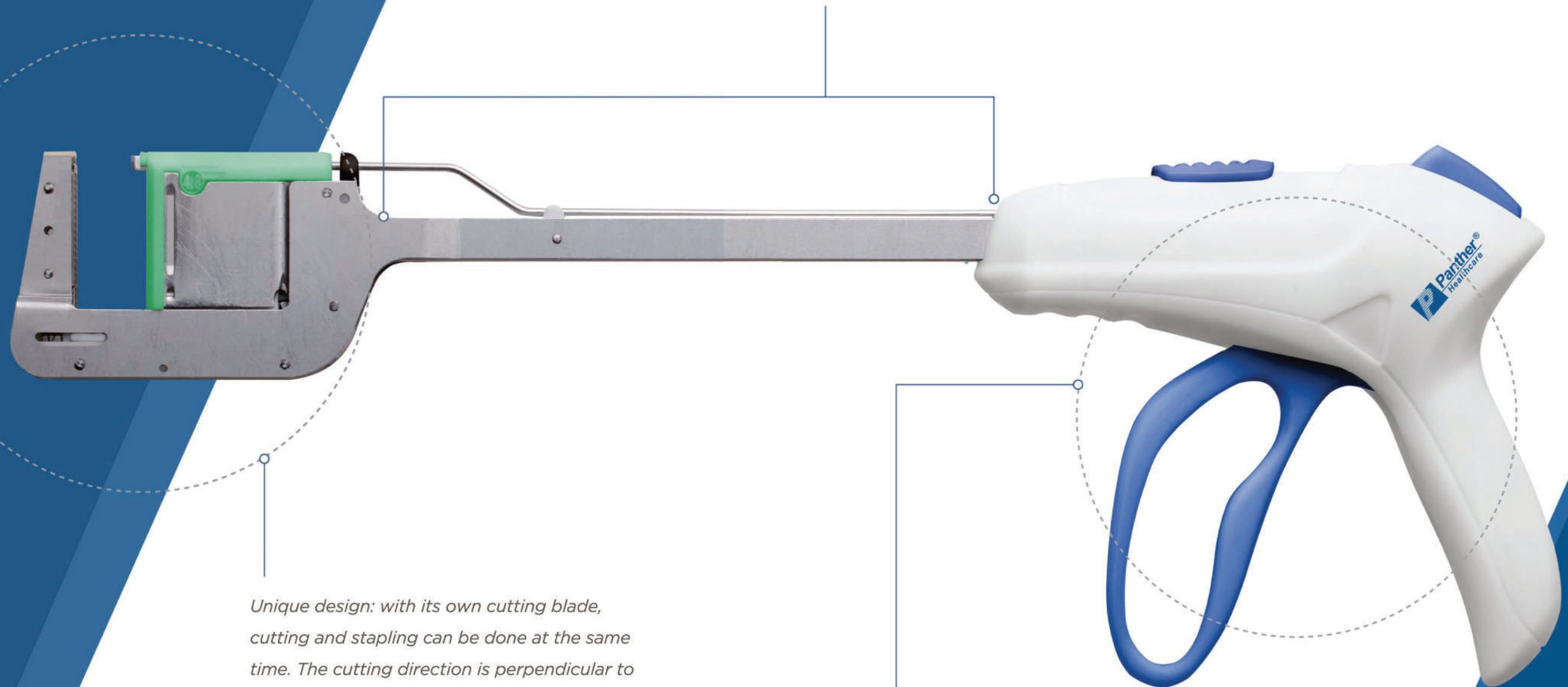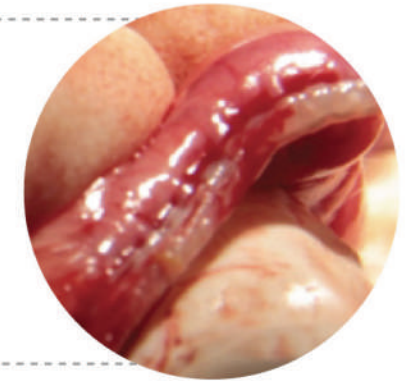

Unique design: with its own cutting blade, cutting and stapling can be done at the same time. The cutting direction is perpendicular to stapling; first stapling then cutting ensures a smoother cut on thicker tissue.

Ergonomic handle design requires less force to fire.